

Operation of Firesting Fiber-Optic Oxygen Meter 2 Channel Software

Date: 2024-09-13
Tags: [O2](#) [AE](#) [Firesting](#) [Software](#)
Category: Protocol
Created by: Alexander Eith

Goal

In situ O2 measurment using the [Equipment - Firesting Fiber-Optic Oxygen Meter 2 Channel](#)

Prerequisites and preparation

Reaction setup

[Equipment - Firesting Fiber-Optic Oxygen Meter 2 Channel](#)

Laptop with pyroscience software

Note: work carefully! Handle espacially the glas fiber cables with care!

Not for use in HTE setup

Steps

Step number	Step description
1	Build the desired reaction setup with either the Equipment - NS14 Schlenk tube for sensor cap or the AE-290: Liquid-phase calibration of Trace Range Robust Oxygen Probe
2	Place the Firesting onto a Kimtech wipe outside the irradiation setup
3	Connect the glas fiber cables / trace robust probe / PT100
4	Start the Laptop and log in (password lab laptop: oxygenevolution)
5	Connect the Firsting with the laptop using an USB cable
6	start: Firesting workbench (?)
7	This window opens: Screenshot-2024-10-09-150309.png
8	To show the graphs: click on graph symbol in the table with channel 1...
9	Check for correct calibration file: Screenshot-2024-10-15-091154.png , click on it, than choose correct calibration file
10	To start logging: press big red button next to Logging not active
11	This window opnes: Screenshot-2024-10-09-150707.png
12	Give appropriate name to sample

13	Choose correct folder (the one in the foto: Screenshot-2024-10-09-150707.png)
14	Press start logging
	To stop measurement:
15	Press big green button, for the appropriate screenshot of the graphs: zoom out completely

Linked experiment

- [AE-290: Liquid-phase calibration of Trace Range Robust Oxygen Probe](#)

Linked resources

Equipment - [NS14 Schlenk tube for sensor cap](#)

Equipment - [Firesting Fiber-Optic Oxygen Meter 2 Channel \(Firesting 1\)](#)

Attached files

Screenshot-2024-10-15-091154.png
sha256: de30fac7010f93b518dac39025037d14b306f0327b94829e1c51f5ef85f48eee



Screenshot-2024-10-09-150436.png
sha256: 12dfe687ba2f4b5a6a211be631fab8f59351f11077a80d220a8f40b2ce1757f6



Screenshot-2024-10-09-150410.png
sha256: c9c1c2ab4143ed7ba45cf6a2845ce50df10ce52f1509e9bb4f1245b4f56aedc9




Screenshot-2024-10-09-150349.png
sha256: 21beae9375605b94e4c543b6cb2065cb9d954ebe2dd698437c44718b58d3bc0d



Screenshot-2024-10-09-150309.png

sha256: 40811810f53b83f06b3c27d15784e185b01b658b71b2bc6d261a7bba661d4f6a

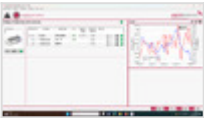

Screenshot-2024-10-09-150707.png

sha256: e7ac662a518c58f2798fbd1470fd246ee2de211f03d5b6e57e797684dfbda1e5




Screenshot-2024-10-09-150524.png

sha256: c311a3b3c7b64c83c9189ecd3a1a57ee5dd632a5aad78e7dc679ba1d262bfb15



Screenshot-2024-10-09-150506.png

sha256: c6dfabf3727931aaecbfdacb8d02df14783d172d8952d47e7b1d7c0f0a5ba25a

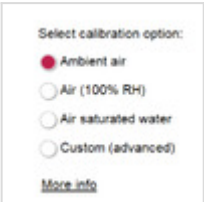

Manual_FS02-C.pdf

sha256: d024469bbe13391701f735ccce7f4402059124bb74f9ef72f1872b477b42290e



Manual_Oxygen_Sensors.pdf

sha256: 1acfa6c06fb26c9f3afbd268dd64319cf17bbae20fcf6389bfb15f38858b0079



Manual_Optical_Temperature_Sensors.pdf

sha256: 5be2676f3c83f0485c5df38bd0c68648104862b4b44147f72f712621fb53aa3b



Manual_Pyro_Workbench-DataInspector.pdf

sha256: 21614a107dcdc1398052674231d23faa15193c75ab8079154c8e14b7cd6d8889



Unique eLabID: 20240913-495002969f505d49021189bc476d408cb6b12fc1

Link: <https://elab.water-splitting.org/database.php?mode=view&id=143>