

NB-321: Modification of EA-354 (SrTiO₃:Al, upscaled batch) with Rh, Cr oxide cocatalyst (0.05 wt%)

Date: 2025-10-27


Tags: NB synthesis Furnace Muffle
Furnace SrTiO₃ RhCl₃ Osterloh
RhxCrO₃ Al:SrTiO₃ RhCrO₃:Al:SrTiO₃
Cr(NO₃)₃

Category: SrTiO₃

Status: Done

Created by: Nadzeya Brezhneva

Reaction scheme/sample structure

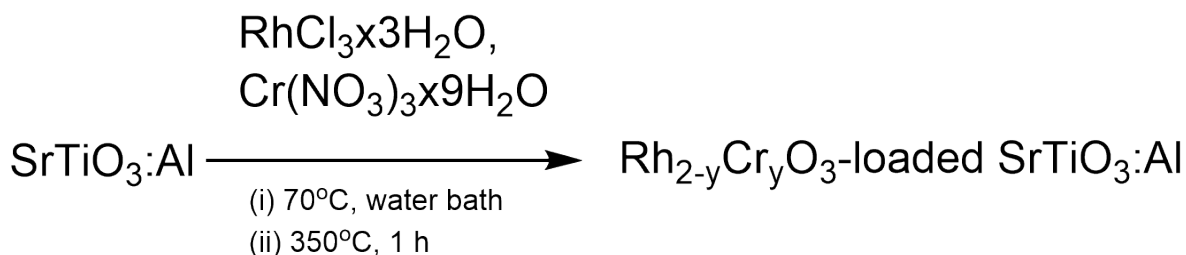

ChemDraw File (linked): [NB-189-RhCr-oxide-loaded-Al-SrTiO3.cdxml](#)

Literature/reference experiments

Literature	https://doi.org/10.1039/C9EE00310J
Reproduction	/
Similar experiments	SrTiO₃ - EA-358: Modification of Al:SrTiO₃ (EA-354) via deposition of Rh, Cr oxide co-catalyst, 350°C, 1h, Upscaling (3.33x) SrTiO₃ - EA-359: Modification of Al:SrTiO₃ (EA-354) via deposition of Rh, Cr oxide co-catalyst, 350°C, 1h

Reagents

Name	CAS Number / Experiment Number	Inventory number	Amount	Equivalents	Mass _{theo} [mg]	Mass _{exp} [mg]	Molar mass [g/mol]	Density [g/ml]	Volume [ul]
Al:SrTiO ₃ (EA-354)	SrTiO₃ - EA-354: Preparation of Al:SrTiO₃ (using EA-352-SrTiO₃-molten-salt and SrCl₂ as flux) at 1000°C, 10 h, Osterloh route, upscaling x15, filtration with PVDF	/	0.81 mmol	1	150	151.03	183.49	/	/
9.87 mM RhCl ₃ solution (EA-357)	#EA-357	/	0.73 umol	0.05 wt% Rh ^[a]	/	75	209.26	/	75
28.8 mM Cr(NO ₃) ₃ solution (EA-357)	#EA-357	/	1.44 umol	0.05 wt% Cr ^[b]	/	50	238.011	/	50
milli-Q water	/	/	0.19 mol	/	/	/	18	0.998	3875 ^[c]

[a] - Rh content in the final product

[b] - Cr content in the final product

[c] - for suspension preparation

Excel sheet for reagent calculation

[NB-321-calculations.xlsx](#)

Furnace Parameters

Equipment - Muffle furnace Nabertherm LT 15/11/P330 (AWZ 304 lab)

Protocol - Muffle furnace Nabertherm GmbH LT 15/11/P330 (Lab AWZ 304)

Temperature/time parameters

Used zone or charge sensor	Zone
Used delayed start	/
Used automatic/manual/extended holdback	Automatic
The temperature band entered for manual/extended holdback (°C)	/
End time [min], relative to start of program	/

Segments

	Target Temperature (°C)	Duration (h)	Rate (°C/h)	Temperature band (°C)	Description of the segment	Observations
First segment	350	0:33	600	/	Increase	/
Second segment	350	1:00	/	/	Hold	/
End segment	0	/	/	/	Natural cooling down	/

Procedure/observations

For transfer of precise liquid amount, Eppendorf pipettes were used (for 10 - 100 μL : 100 μL Eppendorf, above 100 μL : 1000 μL Eppendorf).

Date	Time	Step	Observations	Pictures/Files
27.10.2025	18:40	Weighing EA-354 in a 25 mL Schott glass beaker.	Slightly creamy solid.	/
	ca. 18:45	Placing the glass stirring bar (8 mm) inside the beaker.	/	/
	18:46	Addition of 3875 μL of H_2O ($3 \times 1000 + 875$) to the beaker.	Creamy, slightly pinkish suspension.	20251027_184656-after addition of water.jpg
	ca. 18:48	Placing the beaker with the suspension into the water bath (currently, at room T).	/	/
	18:49	Switching the stirring (400 rpm).	/	/
	18:50	Fast addition of 75 μL of 9.87 mM RhCl_3 solution to a beaker under stirring.	/	/
	18:52	Fast addition of 50 μL of 28.9 mM $\text{Cr}(\text{NO}_3)_3$ solution to a beaker under stirring.	/	/
	18:53	Switching the heating mode on (70 $^\circ\text{C}$, precise mode).	/	20251027_185519-start of heating.jpg
	18:53-20:30	Evaporating of water under constant stirring at 70 $^\circ\text{C}$.	/	/

	20:30	Stop stirring, since all water has been evaporated, removing the stirring bar from the beaker.	In general solid seems quite homogeneous	20251027_203107-after evaporation.jpg
	ca. 20:35	Transfer of the solid to the quartz crucible using Smartspatula, covering the crucible with lid and Al foil before further use.	/	20251027_203754-solid transferred to crucible before calcination.jpg 20251028_004550-crucibles covered with Al foil.jpg
28.10.2025	ca. 10:35	Loading the crucible with the dried sample in a muffle furnace Equipment - Muffle furnace Nabertherm LT 15/11/P330 (AWZ 304 lab) .	/	20251028_103723-crucibles before calcination.jpg 20251028_103806-crucibles inside the furnace.jpg
	10:40	Start of the program (350 °C, 1 h, 10 °/min).	/	20251028_104132-start of calcination.jpg
	19:20	Removal of the sample from the furnace.	T = 163 °C when opening the lid Slightly grey solid.	/
	20:25	Collecting the sample from the beaker, transfer to a 4 mL vial with a screw cap, weighing. m = 148.40 mg	NB-321-0.05 wt% Rh,Cr-oxide loaded Al:SrTiO₃ Slightly grey solid.	20251028_203010-final product.jpg

Product characterization

Sample	Mass [mg]	Purity	Mass _{pure} [mg]	Amount [mol]	Yield [%]	Description	Image	Storage location
NB-321-0.05 wt% Rh,Cr-oxide loaded Al:SrTiO₃	148.40	/	/	/	98.26	Slightly grey solid.	20251028_203010-final product.jpg	Lab E004 (CEEC II), SSC cabinet (SrTiO ₃ -6 container)

Results

Modification of Al:SrTiO₃ sample (**EA-354**, upscaled batch) with Rh, Cr oxide cocatalyst (0.05 wt% Rh and Cr in final product) was performed. The final product (**NB-321-0.05 wt% Rh,Cr-oxide loaded Al:SrTiO₃**) is represented with slightly grey solid, m = 148.40 mg, yield 98.26 %.

Linked resources

Equipment - [Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I, \(Matilda\)](#)

Protocol - [Heat treatment using Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I](#)

Attached files

NB-321-calculations.xlsx

sha256: 950437204a986eb1b69ef05c2814c496b5f0d00f75d0a5fabe98ef62040c1d53

20251028_203010-final product.jpg

sha256: 9e79828cce4643361b3b56bb20ba24b72c643bf0be3f5b6749858b231c5c7d81



20251027_184656-after addition of water.jpg

sha256: 94281b634889230df3cc7c66c0d7d7af6fd7441b30cc58fad1dbdaf10248bd4c



20251027_203754-solid transferred to crucible before calcination.jpg

sha256: 91ebbc36a9407f2818b688ee16c7301496092527b6b94e77b8f3d77126840251




20251027_203107-after evaporation.jpg

sha256: 6d6ce7b4c2f8704a1ee38bb86204f38af3dff0bf2430c89a75aca2ed2393c2ad



20251027_185519-start of heating.jpg

sha256: 2aba07e13bc87192e6de9d8f34d5d67373eb5cefe64d15008588493d1b44487f

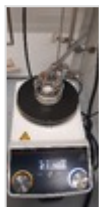

20251028_103723-crucibles before calcination.jpg

sha256: ebcc76b2879f7111363876d90b734899331703ddacdba714dbf9dc6270fa50df



20251028_004550-crucibles covered with Al foil.jpg

sha256: ba59eb481be8c977056f254473977ff960308f6a708c69ce38ecbc2a09f21bb8



20251028_104132-start of calcination.jpg

sha256: e0d01050067f72de0367930d55e595be053c63016e421032e8b04b13483824a7



20251028_103806-crucibles inside the furnace.jpg

sha256: dc978cb64016b945fe874a679a8b1572618ee806ca3c1d5e65e300ab39836e8a



Unique eLabID: 20251027-15a002afdbbccd8d1dca3e503588c4148ce7005b
Link: <https://elab.water-splitting.org/experiments.php?mode=view&id=3263>