

EA-351: Calcination of SrCO₃, 300 C, 1 h, Upscaling (3x)

Date: 2025-10-01

Tags: EA Furnace Muffle Furnace SrTiO₃
SrCO₃ Calcination Upscaling

Category: SrTiO₃


Status: Done

Created by: Ebrahim Abedini

Objective

Calcination of SrCO₃ at 300°C, for 1h to make the material ready for SrTiO₃ synthesis. 3 times upscaling

Reaction scheme/sample structure

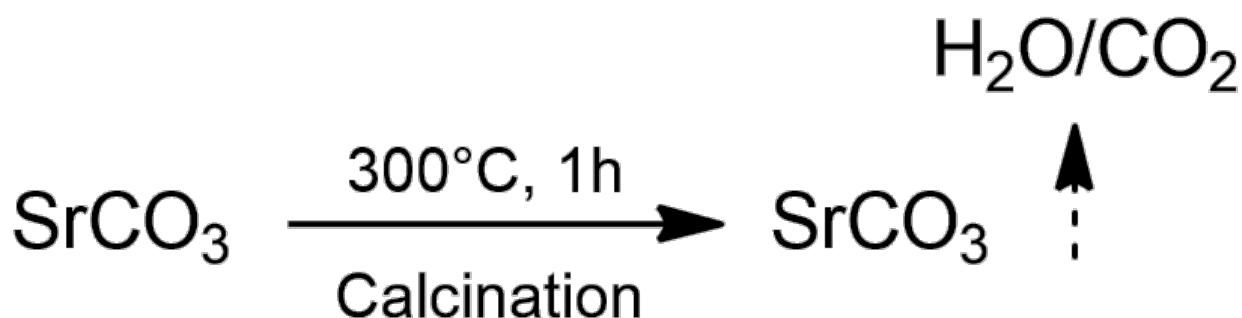

ChemDraw File (linked): [EA-351.cdx](#)

Literature/reference experiments

Literature	https://doi.org/10.1039/C9EE00310J
Reproduction	/
Similar experiments	SrTiO₃ - NB-265: Calcination of SrCO₃, 300 C, 1 h I

Reagents

Name	CAS Number / Experiment Number	Inventory number	Amount [mmol]	Equivalent s	Mass _{theo} [g]	Mass _{exp} [g]	Molar mass [g/mol]	Density (g/ml)	Volume [ml]	Concentration [mM]
Strontium Carbonate (SrCO ₃), TCI	1633-05-2 20250929_1802-SrCO3.jpg	C141117	62.99	/	3 × 3.10 = 9.30	9.31	147.63	3.5	/	/

Furnace Parameters

Equipment - Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I, (Matilda)

Protocol - Heat treatment using Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I

Temperature/time parameters

Used zone or charge sensor	Zone
Used delayed start	/
Used automatic/manual/extended holdback	automatic
The temperature band entered for manual/extended holdback (°C)	/
End time [min], relative to start of program	202

Segments

Program 7	Target Temperature (°C)	Duration (h)	Rate (°C/h)	Temperature band (°C)	Description of the segment	Observations
First segment	300	00:22	600	/	Increase	Initial temperature: 86°C
Second segment	300	01:00	/	/	Holding	/
End segment	/	/	/	/	Natural cooling	End temperature (before opening): 136°C

Procedure/observations

Date	Time	Step	Observations	Pictures
29.09.2025	18:06	Weighing 9.31g SrCO ₃ on a weighing bowl,	White clumps	20250929_1806-weighing the sample.jpg
	18:09	Transferring the sample to an agate mortar	White clumps	20250929_1809-sample before mortaring.jpg
	18:11	The sample was mortared for 2 min.	Broken white clumps, more fine	20250929_1811-sample after mortaring.jpg
	18:13	The mortared sample was transferred into a porcelain crucible ()	/	20250929_1813-the mortared sample before calcination.jpg

	18:21	The crucible with lid was transferred into the furnace Equipment - Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I, (Matilda) . The program was designed according to Protocol - Heat treatment using Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I and saved on P=7.	/	20250929_1821-the crucible inside the furnace.jpg
	18:23	Start of the heating program (EA-351-P=7).	Initial temperature: 86°C	20250929_1823-start of the program.jpg
	21:45	End of the heating program.	End temperature (before opening): 136°C	20250929_2145-End of the program.jpg
	21:49	The crucible was removed from the furnace.	White powder, more puffy	20250929_2149-the sample after calcination.jpg
	21:53	The calcined SrCO_3 was weighed on a weighing bowl after the calcination.	White powder	20250929_2153-weighing the calcined SrCO_3.jpg
	21:58	The weighed sample was transferred into a 10ml snap-cap vial and covered with Al foil. Named: EA-351-Calcined SrCO_3	White powder EA-351-Calcined SrCO_3	20250929_2158-EA-351-Calcined SrCO_3.jpg
		Note: While weighing the sample, the amount for the next step was collected and used for synthesis of SrTiO_3	/	/

Product characterization

Sample	Mass [g]	Purity	Mass _{pure} [g]	Amount [μmol]	Yield [%]	Description	Image	Storage location
EA-351-Calcined SrCO_3	9.24	/	/	/	99.24	White powder	20250929_2158-EA-351-Calcined SrCO_3.jpg	Equipment - Safety storage cabinet - CEEC II Lab E004, shelf A, Ebi container

Results

Calcination of the SrCO_3 was performed at 300°C, for 1h. The yield was 99.24%.

Linked experiment

[SrTiO3 - NB-265: Calcination of \$\text{SrCO}_3\$, 300 C, 1 h I](#)

Linked resources

Equipment - [Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I, \(Matilda\)](#)

Equipment - [Advanced irradiation chamber V1.0 I](#)

Equipment - [Advanced power measurment chamber V1.0 I](#)

Equipment - [Safety storage cabinet - CEEC II Lab E004](#)

Protocol - [Heat treatment using Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I](#)


Attached files

EA-351.cdx

sha256: 6201ae54674d9865fdab4122a0e62d503db2f3cadbc5859983d80a72d7e3f577

EA-351.png

sha256: 88b2d3cb60c6970bcb37a1ff25659fc0d8669153c01a9c73fca79a9e4a5f13a9

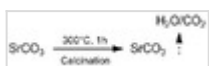

20250929_2158-EA-351-Calcined SrCO3.jpg

sha256: ed3d0f86040345cbf00453e2651be886e067ab4778fcfe11b40481eb8f382d58



20250929_2145-End of the program.jpg

sha256: 6e8ab96832269bcfd1e6531feb67c4fbc11094c86530ce28031e5d010de734db



20250929_2149-the sample after calcination.jpg

sha256: 1284bd142bc76a117ee79f309faa63601a1d2f0d188dbb086141e2b892300389



20250929_2153-weighing the calcined SrCO_3 .jpg

sha256: a41555fe9da8aea8d6e893be3cb8ac824f0c4722391ad3749678f1356cf9882b



20250929_1823-start of the program.jpg

sha256: 8ad7c3a5efa2a8f0aa36c869439239f0ed154efef48aeba1d326cac534aa7953




20250929_1821-the crucible inside the furnace.jpg

sha256: eca125f736e46ad6b28e3df818f7801a68de5338fd16ff4a1a39d2f89d629f65



20250929_1809-sample before mortaring.jpg

sha256: a5de349bd3943fe2f44f3ed3988ff0869a126165d3d2147d844a9861a36a416d

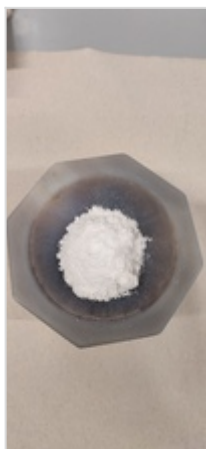

20250929_1811-sample after mortaring.jpg

sha256: 4f4fddd43ab3ed77b7858666dfa09033ffd75851fe8811d812bafc65a433e0a6




20250929_1813-the mortared sample before calcination.jpg

sha256: 7cc884d6f39e632d28133635107ea7fbe9cc9a4723cf3824d3c8ade3ca9c0ff0




20250929_1802-SrCO3.jpg

sha256: 0a7e6f0b2c2e8fc4787a175c83210b381f13587ee3537973a8c40f45995d7710

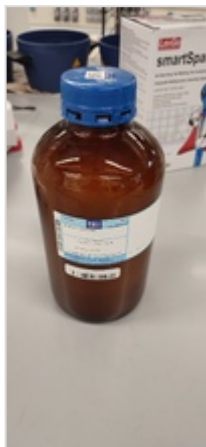

20250929_1806-weighing the sample.jpg

sha256: 95b421ca0205eab69e8f5fe2c3f6b4acd524e6cfd11bc9868ed2646767771ab4

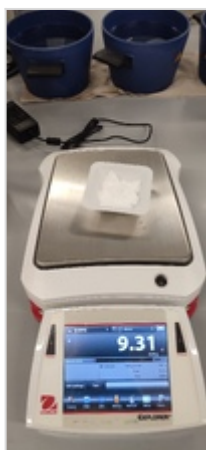

Unique eLabID: 20251001-5091b6ab46e10d50e540c5abf9c24496e5884df9
Link: <https://elab.water-splitting.org/experiments.php?mode=view&id=3075>