

EA-352: Synthesis of SrTiO₃ from SrCO₃ and TiO₂, 1000°C, 10 h, Osterloh route, Upscaling (2x)

Date: 2025-10-01

Tags: Future synthesis EA Furnace
Muffle Furnace SrTiO₃ Osterloh TiO₂
Upscaling

Category: SrTiO₃


Status: Done

Created by: Ebrahim Abedini

Objective

Synthesis of SrTiO₃ from SrCO₃ and TiO₂ (P-25) via solid state approach at 1000°C, 10h. 2 times upscaling.

Reaction scheme/sample structure

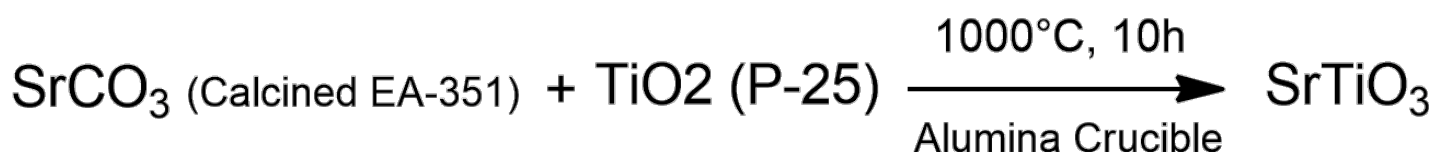

ChemDraw File (linked): [EA-352.cdx](#)

Literature/reference experiments

Literature	https://doi.org/10.1039/C9EE00310J
Reproduction	/
Similar experiments	SrTiO₃ - NB-283: Preparation of SrTiO₃ from SrCO₃ and TiO₂, 1000 C, 10 h, batch V

Reagents

Name	CAS Number / Experiment Number	Inventory number	Amount [mmol]	Equivalents	Mass _{theo} [g]	Mass _{exp} [g]	Molar mass [g/mol]	Density (g/ml)	Volume [ml]	Concentration [mM]
Calcined SrCO ₃ (EA-351)	SrTiO₃ - EA-351: Calcination of SrCO₃, 300 C, 1 h, Upscaling (3x)	/	41.99	1	2 × 3.10 = 6.20	6.21	147.63	3.5	/	/
TiO ₂ , Aeroxide® P25, Thermo Scientific	13463-67-7	C137952	42.06	1	2 × 1.68 = 3.36	3.37	79.87	/	/	/

Work-up and Analytical Reagents

Name	CAS Number / Experiment Number	Inventory number	Mass _{exp} [g]	Volume [ml]	Concentration [M]
------	--------------------------------	------------------	-------------------------	-------------	-------------------

Ethanol	64-17-5	/	/	1	/
---------	---------	---	---	---	---

Furnace Parameters

Equipment - Muffle furnace Nabertherm LT 15/11/P330 (AWZ 304 lab)

Protocol - Muffle furnace Nabertherm GmbH LT 15/11/P330 (Lab AWZ 304)

Temperature/time parameters

Used zone or charge sensor	Zone
Used delayed start	/
Used automatic/manual/extended holdback	automatic
The temperature band entered for manual/extended holdback (°C)	/
End time [min], relative to start of program	861

Segments

Program 2	Target Temperature (°C)	Duration (h)	Rate (°C/h)	Temperature band (°C)	Description of the segment	Observations
First segment	1000	01:38	600	/	Increase	Initial temperature: 20°C
Second segment	1000	10:00	/	/	Holding	/
End segment	/	/	/	/	Natural cooling	End temperature (before opening): 217°C

Procedure/observations

Date	Time	Step	Observations	Pictures
29.09.2025	21:55	Weighing 6.21 g of calcined SrCO ₃ (EA-351) in a weighing bowl.	White clumpy powder	20250929_2155-weighing SrCO3.jpg

	22:03	Weighing 3.37 g of TiO_2 (P-25) in a weighing bowl.	White very fine powder	20250929_2203-weighing TiO2.jpg
	22:06	The weighed samples were transferred into an agate mortar (d: 10cm).	White powders	20250929_2206-transferring the weighed samples into one agate mortar.jpg
	22:12 - 22:20	The weighed SrCO_3 and TiO_2 in the mortar were mixed while mortaring. During the mixing time, 250 μl of distilled EtOH (from the washing bottle) was added to the mixture then got mixed with the powder with a smart spatula. This procedure was repeated 4 times. In total 1ml of EtOH was used.	Addition of EtOH helped the mixture to be more homogenized. After mixing/mortaring the mixture turned to be more fine and less puffy.	20250929_2212-addition of 250 μl EtOH to the mixture-1.jpg 20250929_2214-addition of 250 μl EtOH to the mixture-2.jpg 20250929_2216-addition of 250 μl EtOH to the mixture-3.jpg 20250929_2217-addition of 250 μl EtOH to the mixture-4.jpg 20250929_2220-the mixture before transferring to crucible.jpg
	22:22	The mixture was transferred into a 150ml alumina crucible.	White powder	20250929_2222-transferring the mixture into the alumina crucible.jpg
	22:24	The mixture inside the crucible was pressed with the agate pestle (without applying force and just with help of weight of the pestle).	Pressed and evened white powder	20250929_2224-pressing the mixture in crucible with pestle.jpg
	22:31	The crucible with lid was transferred to AWZ lab 304 Equipment - Muffle furnace Nabertherm LT 15/11/P330 (AWZ 304 lab) . The heating program was designed according to Protocol - Muffle furnace Nabertherm GmbH LT 15/11/P330 (Lab AWZ 304) , and saved on P=2. The heating program (EA-352-P=2) was started.	Initial temperature: 20°C	20250929_2231-crucible inside the furnace.jpg
30.09.2025	00:20	The temperature of the furnace was checked during the heating phase. Note: the furnace can not be heated with 600 (°C/h).	T: 887°C Because of the size of the furnace (15L), the heating takes more time than 10°C·min ⁻¹	20250930_0020-checking the temperature of the program during heating phase.jpg
	12:39	The temperature of the furnace was checked during the cooling phase.	T: 721°C	20250930_1239-checking the temperature of the furnace in cooling down phase.jpg
	23:52	After cooling down the furnace, turned off and opened.	End temperature (before opening): 217°C	20250930_2352-End of heating program.jpg

	00:38 - 00:39	The crucible was opened and the formed SrTiO ₃ was checked.	White big clump with the shape of bottom of the crucible	20251001_0038-formed SrTiO ₃ in crucible after heating program-1.jpg 20251001_0039-formed SrTiO ₃ in crucible after heating program-2.jpg
	00:40	The formed SrTiO ₃ was weighed in a weighing bowl.	White big clump	20251001_0040-Weighing SrTiO ₃ .jpg
	00:42	With the help of a smart spatula, the big clump was broken into smaller pieces.	White soft clumps	20251001_0042-breaking the big clump into smaller pieces.jpg
	00:47	The weighed SrTiO ₃ was collected after breaking the big clump into a 20ml snap-cap vial and covered with Al foil. Named: EA-352-SrTiO₃	White soft clumps EA-352-SrTiO₃	20251001_0047-SrTiO ₃ .jpg
	00:47	The formed SrTiO ₃ of this batch was compared with SrTiO₃ - NB-283: Preparation of SrTiO₃ from SrCO₃ and TiO₂, 1000 C, 10 h, batch V to see the color differences between two batches.	This batch: White NB-283: more creamy/grey	20251001_0047-comparison of this batch with NB-283-1.jpg 20251001_0047-comparison of this batch with NB-283-2.jpg

Product characterization

Sample	Mass [g]	Purity	Mass _{pure} [g]	Amount [μmol]	Yield [%]	Description	Image	Storage location
EA-352-SrTiO ₃	7.64	/	/	/	98.98	White soft clumps	20251001_0047-SrTiO ₃ .jpg	

Results

Synthesis of SrTiO₃ using SrCO₃ and TiO₂ (P-25) was done. The colour of the product was whiter in comparison with previous batches. The mass yield: 79.7%. The theoretical yield: 98.98%.

Future recommendations

Old procedure	Problem	Suggested new procedure
Using Equipment - Muffle furnace Nabertherm LT 15/11/P330 (AWZ 304 lab)	Long time of cooling down	Using Equipment - Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I, (Matilda) instead.

Linked experiment

SrTiO₃ - [EA-351: Calcination of SrCO₃, 300 C, 1 h, Upscaling \(3x\)](#)

Linked resources

Equipment - [Muffelofen L3/11/P320, Nabertherm GmbH, Lab 106, CEEC I, \(Matilda\)](#)

Equipment - [Advanced irradiation chamber V1.0 I](#)

Equipment - [Advanced power measurment chamber V1.0 I](#)

Equipment - [Muffle furnace Nabertherm LT 15/11/P330 \(AWZ 304 lab\)](#)

Equipment - [Manual irradiation setup](#)

Protocol - [Muffle furnace Nabertherm GmbH LT 15/11/P330 \(Lab AWZ 304\)](#)

Attached files

EA-352.cdx

sha256: 89aea3de229e2f4eca354bc3f5e19188f9d34c2cf0d06d835e74ce8cac1379e3

EA-352.png

sha256: 37c623ae90cc4a81ffcaa54744deaf1113d664e53b315a91cc830df0f0140ad5



20250929_2155-weighing SrCO₃.jpg

sha256: 4bb36ff69f36c8c7b357de4b6dbb5dce13eec853bec93a9804097b011ee8cdc6



20250929_2203-weighing TiO₂.jpg

sha256: fa5f65b399dc0464551c0667b0d10fbd2a17a2439962bee7ba67bd863b3ca6b7



20250929_2206-transferring the weighed samples into one agate mortar.jpg

sha256: 139f28e9fe3f50adaff256e8c39ff754775d21c7593a4865eb3f6930a3bcb5d6




20250929_2212-addition of 250µl EtOH to the mixture-1.jpg

sha256: ce4af6332b2e1c34d6a72830bc5b988cd4cdc7811c498305392d43ed912779f6




20250929_2214-addition of 250µl EtOH to the mixture-2.jpg

sha256: 019d4dafb56e46cc32b6c5c1a66f97f84e869d77b454f0f590d2621cb67fcb9c

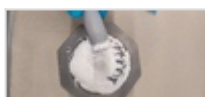

20250929_2216-addition of 250µl EtOH to the mixture-3.jpg

sha256: 9ce9e8a40762e4b50b5f90f6cfa55c7e5b0772b861be982194ae76c638f66119

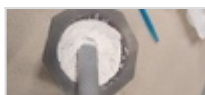

20250929_2217-addition of 250µl EtOH to the mixture-4.jpg

sha256: 680f480c26efb1b9e5ab918d73dd284f8db28f4c02a3d22a9edae4dc02d68c2c



20250929_2220-the mixture before transferring to crucible.jpg

sha256: 9b4067920c79ec185642faff9738924834d7a5e6e4875513159b0ed6d864c60f




20250929_2224-pressing the mixture in crucible with pestle.jpg

sha256: f060e82ff927b895445f15fce034e0be67f3133197ccc5c4bf26aab423c10850



20250929_2222-transferring the mixture into the alumina crucible.jpg

sha256: 72aabfe8447f894ed9626ffb0983852fabeeb3347da433ac7a15866f3acbcd6

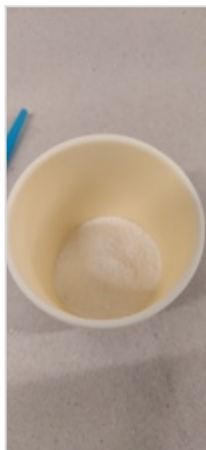

20250929_2231-crucible inside the furnace.jpg

sha256: 0fc26ee819b0354c04ccc9d7eae0f7cd6d9b3d6d9c1702e303a4970388e8deaa



20250930_0020-checking the temperature of the program during heating phase.jpg
sha256: 15132f68052ac28aa6442474d8cf7782e8778f10890a268ac5aafc4ce1583a2



20250930_1239-checking the temperature of the furnace in cooling down phase.jpg
sha256: 82449a56b02acc919c04be353f88096d35d1a22f50fac27c0605dee6dcd75033



20250930_2352-End of heating program.jpg
sha256: 89bcbcd25fb09f6a3b6b8231e9895a3818aba54bb63c8d4ff28fe7b37db133257



20251001_0038-formed SrTiO₃ in crucible after heating program-1.jpg
sha256: 62623bd1b704d651244789e2d64ddfad3ec336f9c5505a4331f7974b606cc9a1



20251001_0039-formed SrTiO₃ in crucible after heating program-2.jpg
sha256: 295aac063d231fa1c06880dcdbf1e9d137054a3aad195161987218fc6bd7ed2b



20251001_0040-Weighing SrTiO₃.jpg

sha256: 392edf03163fbce43f4751d73f08dd89a6c90efc195ef2ae5e68fe041b7d0672




20251001_0042-breaking the big clump into smaller pieces.jpg

sha256: 7164053d51fe4be8eabf2602831d4ea4486dd2a52d3ea6b516220e9457b9a6ec




20251001_0047-comparison of this batch with NB-283-2.jpg

sha256: 715bb9813c48a1a0f1177fb868de6c7c13b9acbe301575a9d0d35ff00641c9a9

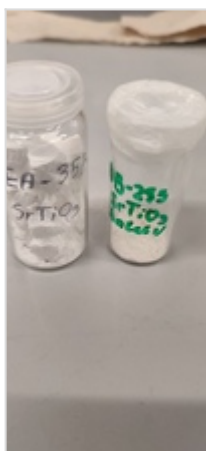


20251001_0047-comparison of this batch with NB-283-1.jpg

sha256: 31cd09f71830b5f0d041d64588c7b30cd93aee23c3cf5852ae72a37654029759

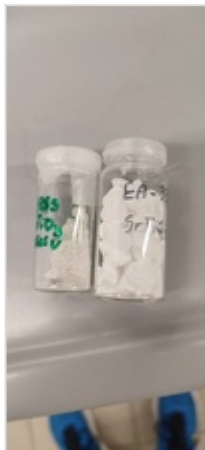

20251001_0047-SrTiO3-molten-salt.jpg

sha256: 85e30ac3211cabe54befb3f2fe0c71b66f0bf1e27520f0127b1040f74af3cf49

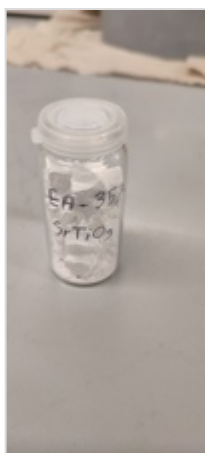

Unique eLabID: 20251001-729a0f21b29fa0b278b770c44e4183869a1e3f69
Link: <https://elab.water-splitting.org/experiments.php?mode=view&id=3074>