

NB-328: Liquid phase H₂ and O₂ of RhCrO_x,Al:SrTiO₃ (NB-323, 0.2 wt% Rh, Cr, 0.5 mg/mL), 365 nm, 50 mW/cm², 20 °C

Date: 2025-10-28

Tags: O₂ Test Calibration Future NB
Firesting O₂ sensor H₂ SrTiO₃ Unisense
RhCrO₃:Al:SrTiO₃ H₂ Sensor
temperature In situ Trace range robust
oxygen sensor photocatalysis Unisense
normal range

Category: SrTiO₃


Status: Need to be redone

Created by: Nadzeya Brezhneva

Objectives

Simultaneous detection of H₂ and O₂ evolution in liquid phase for irradiated suspension of Rh_xCr_yO₃:Al:SrTiO₃ suspension (NB-323, **0.2 wt% Rh,Cr**, 0.5 mg/mL), 365 nm LED, 50 mW/cm², 20 °C.

Reaction scheme

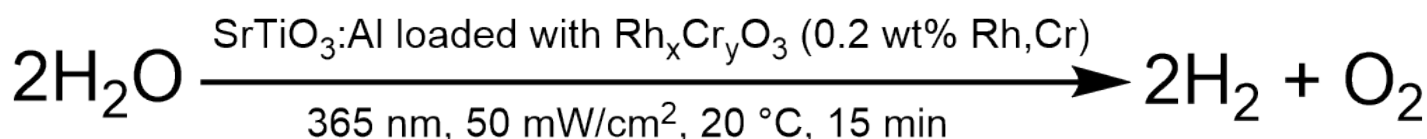

ChemDraw file linked: [NB-328-SrTiO3-photocatalytic H2O splitting.cdxml](#)

Literature/reference experiments

Literature	/
Reproduction	/
Similar experiments	SrTiO3 - NB-316: Liquid phase H₂ and O₂ of RhCrO_x,Al:SrTiO₃ (EA-358, 0.5 mg/mL), 365 nm, 50 mW/cm², 20 °C SrTiO3 - NB-327: Liquid phase H₂ and O₂ of RhCrO_x,Al:SrTiO₃ (NB-321, 0.05 wt% Rh, Cr, 0.5 mg/mL), 365 nm, 50 mW/cm², 20 °C

Reagents

Name	CAS Number / Experiment Number	Inventory number	Amount [mmol]	Equivalents	Mass _{theo} [mg]	Mass _{exp} [mg]	Molar mass [g/mol]	Density (g/ml)	Volume [ml]	Concentration [mM]
milli-Q H ₂ O	/	/	/	/	/	/	/	0.998	25	/
Al:SrTiO ₃ RhCrO _x (NB-323)	SrTiO₃ - NB-323: Modification of EA-354 (SrTiO₃:Al, upscaled batch) with Rh, Cr oxide cocatalyst (0.2 wt%)	/	/	/	12.50	12.54	/	/	/	/

Excel sheet for reagent calculation

/

Irradiation Parameters

Power measurement was performed using [Power Meter - 843-R-USB + 919P-020-12](#) in [Equipment - Advanced power measurment setup V1.0 I](#)

Power measurement was performed in experiment [Prep work - NB-314: Measuring power output of UHP-365 nm #4 with 18A-4 in advanced irradiation setup](#)

	Name
Used Set-up	Equipment - Advanced irradiation setup V1.0 I
Irradiation setup number	Equipment - Irradiation setup 4 (CEEC II, E002)

	Light Source Name	Power Source Name	Wavelength [nm]	Power Setting [mW]	Analog Setting [0.00 - 10.00]
First light source	Light Source - UHP LED 365 nm-4	Power Sources - BLS-18000-1 4	365	56	0.19

Used beam combiner [Name or None]	/
Irradiation distance [cm]	6.5
Thermostat temperature [°C]	20
Stirring speed [rpm]	500
Irradiation start: 1. Firesting [relative to start log] 2. Unisense	1. 601 s 2. 13:13:22
Irradiation stop: 1. Firesting [relative to start log] 2. Unisense	1. 1516 s 2. 13:28:38

O₂/H₂ sensor equipment

	Equipment	Used protocol
Used Firesting	Equipment - Firesting Fiber-Optic Oxygen Meter 2 Channel (Firesting 2)	Protocol - Operation of Firesting Fiber-Optic Oxygen Meter 2 Channel Software
Used O ₂ sensor	Equipment - Robust probe for liquid O ₂ measurment	Protocol - In-situ hydrogen and oxygen measurment in H ₂ /O ₂ reactor
Used H ₂ sensor	Equipment - H ₂ UniAmp Sensor - Normal range - 2.1 x 80 mm needle	Protocol - In-situ hydrogen and oxygen measurment in H ₂ /O ₂ reactor

Procedure/observations

Date	Time	Step	Observations	Pictures/Files
		Calibration from SrTiO ₃ - NB-326: Liquid phase H ₂ and O ₂ of RhCrOx,Al:SrTiO ₃ (EA-358, 0.5 mg/mL), 365 nm, 20 mW/cm ² , 20 °C (reproduction NB-325) was used.	/	/
28.10.2025		Sample preparation		
	20:45	Weighing photocatalyst in a 50 mL vial (the vial was covered with Al foil before further use).	Slightly grey solid	/
29.10.2025	11:30	Addition of 25 mL H ₂ O to the vial via graduated cylinder.	/	/
	11:34-37	The suspension was vortexed for 3 min (Equipment - VWR® VV3, Vortex Mixer, stage 4/6), covered with Al foil before further use.	Seems less optically dense than for the suspension with 0.1 wt% Rh,Cr catalyst EA-358	20251029_113755-suspension after vortex.jpg
		Continue in Protocol - In-situ hydrogen and oxygen measurment in H ₂ /O ₂ reactor from step 6		
	11:40	The suspension was transferred to the reactor using glass pipette (preliminary the vial was manually shaken ca. 15 s) .	/	/
	ca. 11:45	Assembling the setup.	Currently, stopper instead of H ₂ sensor, PT100, PT1000 and O ₂ robust probe are inside the reactor immersed in the liquid phase	20251029_115502-suspension inside the reactor before degassing.jpg

	11:56	Start of O2 logging.	NB-328-Ch2-1	2025-10-29_115632_NB-328-Ch2-1.txt 2025-10-29_115632_NB-328-Ch2-1.png
	11:59	The degassing was started	/	20251029_120001-during degassing of the suspension.jpg
	12:51	Cannula was transferred to gas phase, above the suspension.	/	/
	12:55	H ₂ sensor was added in Ar counterflow.	/	/
	ca. 12:59	The degassing was stopped by removing the cannula and closing the valve.	/	/
	13:01	Stop of O2 logging.	/	/
	13:03	Start of O2 logging.	NB-328-Ch2-2	2025-10-29_130322_NB-328-Ch2-2.txt 2025-10-29_130322_NB-328-Ch2-2.png
	13:03	Start of H2 logging.	NB-328-Logger1	NB-328.ulog NB-328-Logger1.csv NB-328-Logger1.bmp
	13:03-13	Equilibration time.	/	/
	13:13	The irradiation was started	/	20251029_131433-after start of irradiation.jpg
	13:28	The irradiation was stopped.	/	/
	13:28-35	Equilibration time.	/	/
	13:35	Stop of O2 and H2 logging.	/	/

	13:40	Deassembling the setup, cleaning the reactor.	Tips of the sensors and reactor were covered with attached photocatalyst particles. Tip: After preliminary cleaning with sticks, wipes, the residual particles attached to the walls of the reactor could be removed by sonication - fill the reactor with water and place it in ultrasonic bath for ca. 20 s (Eco mode).	/
--	-------	---	--	---

Analysis

Used calibration for Firesting: [20250611-BOLA-fitting-liquid-phase-trace-oxygen-sensor-H2-O2 reactor.ini](#)

Used calibration for UniSense: NB-326-Logger2

Date	Time	Sample name	Analysis method	Analytical device	Solvent	Raw Data	Python script	Processed Data	Comparative Data	Interpretation
29.10.2025	13:03	NB-328-Logger1	electrochemical H2 detection	Equipment - H2 UniAmp Sensor - Normal range - 2.1 x 80 mm needle	water	NB-328.ulong NB-328-Logger1.csv	NB-328-O2 nad H2 curve.py	NB-328-Logger1.bmp NB-328-O2 and H2 curves.png	/	Clean response signal, but sharp increase in the signal after strat and stop of irradiation. H2 value 33.5 uM H2
	11:56	NB-328-Ch2-1	Optical O2 detection	Equipment - Firesting Fiber-Optic Oxygen Meter 2 Channel	water	2025-10-29_115632_NB-328-Ch2-1.txt	/	2025-10-29_115632_NB-328-Ch2-1.png	/	Degassing of the suspension.
	13:03	NB-328-Ch2-2	Optical O2 detection	Equipment - Firesting Fiber-Optic Oxygen Meter 2 Channel	water	2025-10-29_130322_NB-328-Ch2-2.txt	NB-328-O2 nad H2 curve.py	2025-10-29_130322_NB-328-Ch2-2.png NB-328-O2 and H2 curves.png	/	A bit sharp increase in O2 value after start and stop of irradiation, O2 value 27.5 uM.

Results

Simultaneous H₂ and O₂ measurements of irradiated suspension of NB-323 (0.5 mg/mL) in O₂/H₂ photoreactor under 365 nm irradiation (50 mW/cm², 20 °C, 15 min) were performed.

H₂ level at the end of irradiation - 33.5 umol/L, O₂ level - 27.5 umol/L.

Future recommendations

Old procedure	Problem	Suggested new procedure
---------------	---------	-------------------------

/	Sharp increases in O2 and H2 signals - maybe related either with the position of the sensors in the reactor or the optical density.	<ul style="list-style-type: none"> • Adjust the position of the sensors relatively to the light beam. • If necessary, destroy the lumps of the solid photocatalyst before adding water and vortexing.
---	---	---

Linked experiments

SrTiO3 - NB-316: Liquid phase H2 and O2 of RhCrOx,Al:SrTiO3 (EA-358, 0.5 mg/mL), 365 nm, 50 mW/cm2, 20 °C

SrTiO3 - NB-318: Liquid phase H2 and O2 of RhCrOx,Al:SrTiO3 (EA-358, 0.5 mg/mL), 365 nm, 100 mW/cm2, 20 °C

SrTiO3 - NB-323: Modification of EA-354 (SrTiO3:Al, upscaled batch) with Rh, Cr oxide cocatalyst (0.2 wt%)

SrTiO3 - NB-325: Liquid phase H2 and O2 of RhCrOx,Al:SrTiO3 (EA-358, 0.5 mg/mL), 365 nm, 20 mW/cm2, 20 °C

SrTiO3 - NB-326: Liquid phase H2 and O2 of RhCrOx,Al:SrTiO3 (EA-358, 0.5 mg/mL), 365 nm, 20 mW/cm2, 20 °C (reproduction NB-325)

SrTiO3 - NB-327: Liquid phase H2 and O2 of RhCrOx,Al:SrTiO3 (NB-321, 0.05 wt% Rh, Cr, 0.5 mg/mL), 365 nm, 50 mW/cm2, 20 °C

Linked resources

Equipment - [VWR® VV3, Vortex Mixer](#)

Equipment - [Firesting Fiber-Optic Oxygen Meter 2 Channel \(Firesting 2\)](#)

Equipment - [Irradiation setup 4 \(CEEC II, E002\)](#)

Protocol - [Hydrogen Measurement Using GC](#)

Protocol - [Getting hydrogen from hydrogen bottle in CEEC II E014](#)

Protocol - [In-situ hydrogen and oxygen measurment in H2/O2 reactor](#)


Attached files

NB-328-O2 nad H2 curve.py

sha256: 1e5c0cc028f3ae47eec3942074b17fe23d6334b83ed9151d1f6a1ed2e0fa8a12

NB-328-O2 and H2 curves.png

sha256: 83f834eea93fbfc96e1f0cc7f87c00fc8a6556e47bee4a3692c2195412f98605

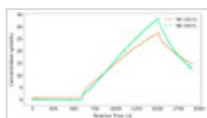

NB-328-SrTiO3-photocatalytic H2O splitting.cdxml

sha256: e3f6550ad862b26ddb6fc1f4ef4ad271ffb39900dfc9788af31e12d3a3e071f2

NB-328-SrTiO3-photocatalytic H2O splitting.png

sha256: 18406fcbf56301b39c7bb7b1d4a3168a1004fa9a9bc42b28670cb24549ef14d1



20251029_115502-suspension inside the reactor before degassing.jpg

sha256: f55b5abac9197b86a9d18f7fccb18821c2acb4f6f720ab38af3078594f0f8ee1



20251029_113755-suspension after vortex.jpg

sha256: 71519753768dbc1c132eec4c229ff561434b34d40d709df8fa2fa5bf0e9e8f9b



20251029_120001-during degassing of the suspension.jpg

sha256: 2d80456974ea2f9f398dff1e62561c6c6f488aa0f087dd9d1a512e84758534fc



20251029_131433-after start of irradiation.jpg

sha256: 3f3aba715bbf4acc59e10c8eae8306bbfaf4c29ba486625f98d94ea295d2e498



2025-10-29_130322_NB-328-Ch2-2.png

sha256: e4d9d9ebf26dbcf0660779751471aa3bbb74ff378eaa0fc92ffd4c988ef62348



2025-10-29_130322_NB-328-Ch2-2.txt

sha256: 9dac6d1589d31db65072964955cfd6cbfaa28bba09771611495bd32669e01a5

2025-10-29_115632_NB-328-Ch2-1.txt

sha256: edacc0f5c4adcbd5a32f831febd56a24be570fc41cf96c16e088e6ffb7156fa5

2025-10-29_115632_NB-328-Ch2-1.png

sha256: dabbeebde233d1fc212f3456f5500d0e165d2ecfbb47b9feb0f15c92201093ca



NB-328-Logger1.bmp

sha256: 7e79f6f42bc8e484a40e4ee5839b34ddeb301fbc54366db6b0af531569e7c7de



NB-328-Logger1.csv

sha256: fa24a3ff3d52143b76d0a512f51efd03c9655939aa73b22b085d46935c5a6fdd

NB-328.ulog

sha256: f688bc0ad3f07c8f6bc6af3c98799fae3f657ac9b38236cd1f8511069c36b8e1



Unique eLabID: 20251028-371c342f1011f30bb49ec0cbb66282c2632d2760
Link: <https://elab.water-splitting.org/experiments.php?mode=view&id=3290>